

Cordless Endomotor

E-CONNECT



5

Reciprocation Modes (M0)

Forward Reciprocation-Four Modes Clockwise&counter clockwise movement



Fwd 150°-Rev 30°|Fwd 180°-Rev 30°

Fwd 210°-Rev 30°|Fwd 250°-Rev 30°

For all continuous rotation files

Reverse Reciprocation-One Mode Counter clockwise&clockwise movement



Fwd 30°-Rev 150°

For reciprocating files



One Piece Contra Angle



Autoclavable



Detachable



Miniature Head



360° Adjustable Contra Angle



Maximum Speed 1000 rpm

RPM

360° Clockwise Rotational Movement

360°

Editable Speed Setting

Editable Torque setting



Forward Direction of Rotation

Reverse Direction of Rotation


Torque Auto Reverse

Editable Memory Programs (M1 to M9)



E-CONNECT × E-PEX

Integrate Endomotor with Apex Locator(Connect Function)

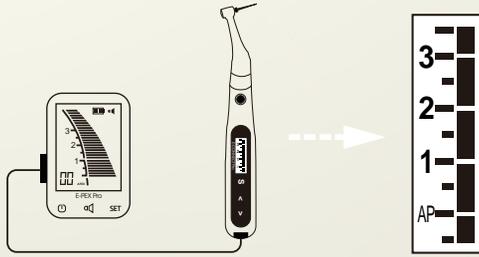

Auto Start & Stop



Apical Reverse



Apical Slow Down



*User can rotate the Endomotor Screen towards left or right

*In connect function Endomotor mirrors Apex Locator's screen



Torque Calibration



Digital OLED Display



Low Noise Motor



Ambidextrous Display



Powerful 1500 mAh Lithium Battery



Integration Option



Screen Cast / Mirroring

Unique Features of E-CONNECT Endomotor