

E-value

Brushless EndoMotor



E-value

Brushless EndoMotor

Contra Angle

Advanced 1:1 contra-angle
One piece design



New contra-angle head



Normal head & neck



360° Rotatable



Brushless Motor

Reliable, Durable, Efficient



Built-in File System

Switch file systems quickly



OLED Color Screen

More comfortable visual experience



Long Battery Life

Heavy capacity
1900mAh lithium battery



E-value

Brushless EndoMotor

Torque Calibration

For precise torque output

Ambidextrous Display

Screen can be rotated, convenience for lefties

Seamless Shell

Easier to clean

Multiple Customizable Modes

Built-in 10 modes that could be set and choose according to actual needs

Built-in File System

① Built-in File System

My Mode	
M1	350rpm
Fwd	2.0N·cm

② Long Press S to enter the Mode Setting interface

Fwd 300rpm 3.0N·cm
My Mode
E-FLEX BLUE
E-FLEX S
E-FLEX

③ Choose the Endo File being used and press S

E-FLEX BLUE	
M1	300rpm
Fwd	3.0N·cm

Note: the doctor could adjust the speed, torque limit and reciprocating angle according to the needs of the case.

Working Modes:


- **Fwd** - Forward Rotation
- **Rev** - Reverse Rotation
- **REC** - Reciprocation Mode

Doctor could set the reciprocating angle and speed level according to actual needs.

- Speed Level range from Level 1 to Level 5
- Angles range from 30° to 370°

- **ATC** - Adaptive Torque Control Mode

After reaching the preset torque limit, motor will reciprocate according to the set angle to improve the efficiency of root canal enlargement.

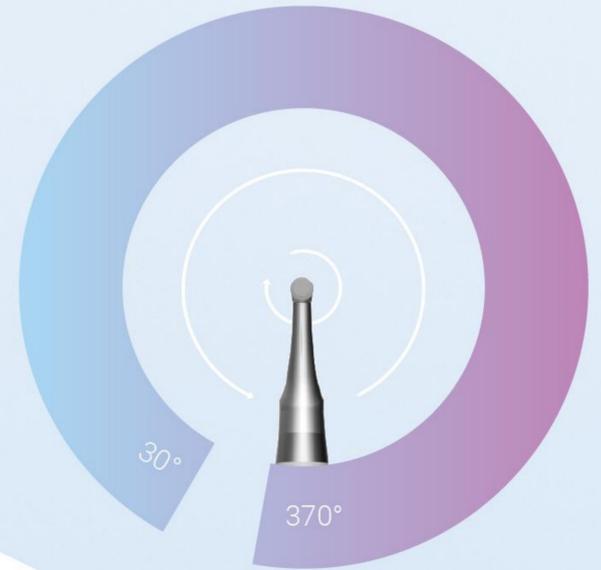

Specifications:

Speed Range	: 120rpm - 1000rpm
Torque Range	: 0.5 - 4N.cm
Memory Programs	: 10 (M1 to M10)
Battery Capacity	: 1900 mAh

